



PROSTHETIC PRODUCTS CATALOGUE


INDEX

HOW TO DETERMINE THE SIZE	3
TRANSTIBIAL SILICONE LINERS.....	4
TRANSTIBIAL TPE Gel LINERS.....	16
TRANSFEMORAL SILICONE LINERS.....	26
KNEE SLEEVE.....	32
CARBON FIBER FEET.....	36
ADAPTERS.....	61

HOW TO DETERMINE THE SIZE


To determine the correct size of the liner:

- Measure the circumference of the residual limb end 4 cm above the distal end of the residual limb.
- To determine the correct size reduce circumference by 0-1.

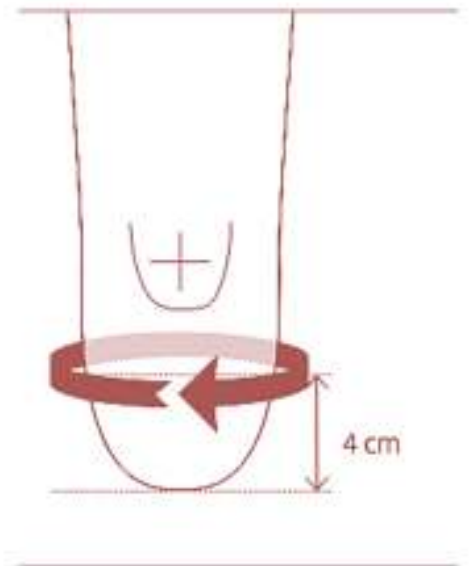

To determine the correct size of the Sleeve:

- Measure the thigh circumference 15 cm above the patella and choose the appropriate size of suspension sleeve.

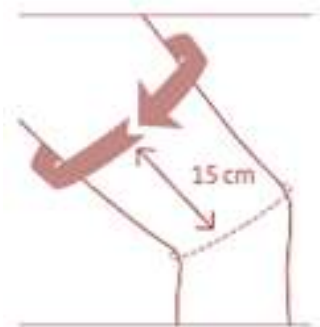

TruSil Pro Liner – Locking



The TruSil Pro liner is designed with very soft silicone gel and with a durable fabric cover. The TruSil Pro liner includes many skin care ingredients to protect the skin and to ensure to the user the maximum comfort.

TruSil Pro liner is indicated for users with a diabetic, dry and allergic skin and with vascular problems.

Features and Benefits:

- Very soft silicone gel to offer the best comfort
- Skin care formulation for care for the skin
- Durable fabric cover for easy doning and doffing
- Flexible distal cup
- Silky internal surface

Indications

Mobility level	Low to High activity
Amputation level	Transtibial

General information

Profile	3 mm – 6 mm
Suspension Method	Locking
Matrix	10 cm
Sizes	22,23.5,25,26.5,28,30,32,34,36,38,40,
Article Number	• FD143xx (3mm profile) • FD146xx (6mm profile)




TruSil Pro Liner – Ordering Guidelines

PROFILE

6 mm uniform

3 mm uniform

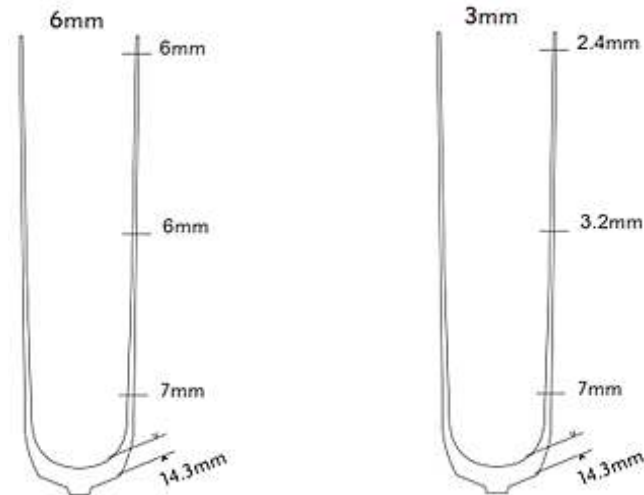

SIZING GUIDE

Measure the circumference of the stump at 4 cm above the distal end.

Measured sized (cm)	3mm Uniform	6mm Uniform
22-23.5 cm	FD14322	FD14622
23.5-25 cm	FD14323.5	FD14623.5
25-26.5 cm	FD14325	FD14625
26.5-28 cm	FD14326	FD14626.5
28-30 cm	FD14328	FD14628
30-32 cm	FD14330	FD14630
32-34 cm	FD14332	FD14632
34-36 cm	FD14334	FD14634
36-38 cm	FD14336	FD14636
38-40 cm	FD14338	FD14638

TruSil Pro Liner – Cushion



The TruSil Pro liner is designed with very soft silicone gel and with a durable fabric cover. The TruSil Pro liner includes many skin care ingredients to protect the skin and to ensure to the user the maximum comfort.

TruSil Pro liner is indicated for users with a diabetic, dry and allergic skin and with vascular problems.

Features and Benefits:

- Very soft silicone gel to offer the best comfort
- Skin care formulation for care for the skin
- Durable fabric cover for easy doning and doffing
- Flexible distal cup
- Silky internal surface

Indications

Mobility level	Low to High activity
Amputation level	Transtibial

General information

Profile	3 mm, 6mm
Suspension Method	Cushion
Matrix	Without
Sizes	22,23.5,25,26.5,28,30,32,34,36,38,40,
Article Number	• FD041xx (3mm profile) • FD042XX (6mm profile)




TruSil Pro Liner – Ordering Guidelines

PROFILE

6 mm uniform

3 mm uniform

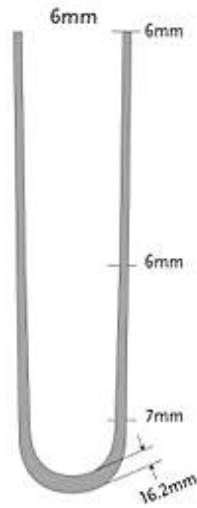

SIZING GUIDE

Measure the circumference of the stump at 4 cm above the distal end.

Measured sized (cm)	3mm Uniform	6mm Uniform
22-23.5 cm	FD04322	FD04622
23.5-25 cm	FD04323.5	FD04623.5
25-26.5 cm	FD04325	FD04625
26.5-28 cm	FD04326	FD04626.5
28-30 cm	FD04328	FD04628
30-32 cm	FD04330	FD04630
32-34 cm	FD04332	FD04632
34-36 cm	FD04334	FD04634
36-38 cm	FD04336	FD04636
38-40 cm	FD04338	FD04638

TruSil Liner – Locking



The TruSil liner is designed with soft silicone gel and with a durable fabric cover. The TruSil liner is indicated for users with a bony residual limbs tenderness to pressure.

Features and Benefits:

- Soft silicone gel to offer the best comfort
- Durable fabric cover for easy doning and doffing
- Flexible distal cup
- Silky internal surface

Indications

Mobility level	Low to High activity
Amputation level	Transtibial

General information

Profile	3 mm, 6mm
Suspension Method	Locking
Matrix	10 cm
Sizes	22,23.5,25,26.5,28,30,32,34,36,38,40,
Article Number	• FD133xx (3mm profile) • FD136xx (6mm profile)




TruSil Liner – Ordering Guidelines

PROFILE

6 mm uniform

3 mm uniform

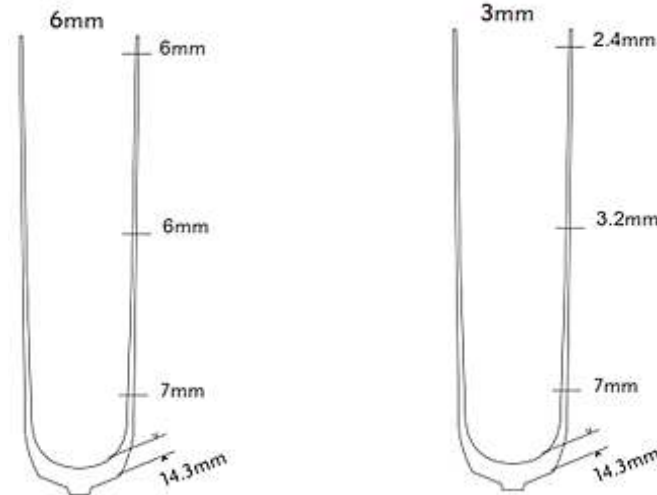

SIZING GUIDE

Measure the circumference of the stump at 4 cm above the distal end.

Measured sized (cm)	3mm Uniform	6mm Uniform
22-23.5 cm	FD13322	FD13622
23.5-25 cm	FD13323.5	FD13623.5
25-26.5 cm	FD13325	FD13625
26.5-28 cm	FD13326	FD13626.5
28-30 cm	FD13328	FD13628
30-32 cm	FD13330	FD13630
32-34 cm	FD13332	FD13632
34-36 cm	FD13334	FD13634
36-38 cm	FD13336	FD13636
38-40 cm	FD13338	FD13638

TruSil Liner – Cushion



The TruSil liner is designed with soft silicone gel and with a durable fabric cover. The TruSil liner is indicated for users with a bony residual limbs tenderness to pressure.

Features and Benefits:

- Soft silicone gel to offer the best comfort
- Durable fabric cover for easy doning and doffing
- Flexible distal cup
- Silky internal surface



Silicone Liners - Transtibial

Indications

Mobility level	Low to High activity
Amputation level	Transtibial

General information



Profile	3 mm, 6mm
Suspension Method	Cushion
Matrix	Without
Sizes	22,23.5,25,26.5,28,30,32,34,36,38,40,
Article Number	• FD033xx (3mm profile) • FD036xx (6mm profile)

TruSil Liner – Ordering Guidelines

PROFILE

6 mm uniform

3 mm uniform

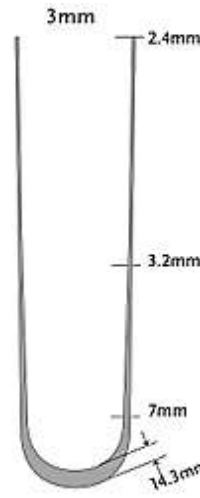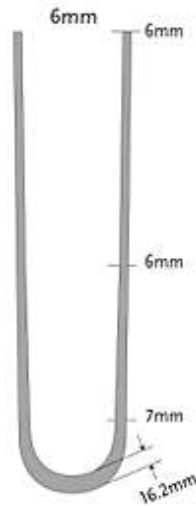

SIZING GUIDE

Measure the circumference of the stump at 4 cm above the distal end.

Measured sized (cm)	3mm Uniform	6mm Uniform
22-23.5 cm	FD03322	FD03622
23.5-25 cm	FD03323.5	FD03623.5
25-26.5 cm	FD03325	FD03625
26.5-28 cm	FD03326	FD03626.5
28-30 cm	FD03328	FD03628
30-32 cm	FD03330	FD03630
32-34 cm	FD03332	FD03632
34-36 cm	FD03334	FD03634
36-38 cm	FD03336	FD03636
38-40 cm	FD03338	FD03638

TruSil Prim Liner – Locking



The TruSil prim liner is designed with firm silicone and with a durable fabric cover. The TruSil Prim liner is for day to day use.

Features and Benefits:

- Firm silicone to offer the best comfort
- Durable fabric cover for easy doning and doffing
- Flexible distal cup
- Silky internal surface



Silicone Liners - Transtibial

Indications

Mobility level	Low to High activity
Amputation level	Transtibial


General information

Profile	3 mm,
Suspension Method	Locking
Matrix	10 cm
Sizes	22,23.5,25,26.5,28,30,32,34,36,38,40,
Article Number	• FD12xx

TruSil Prim Liner – Ordering Guidelines

PROFILE

3 mm uniform

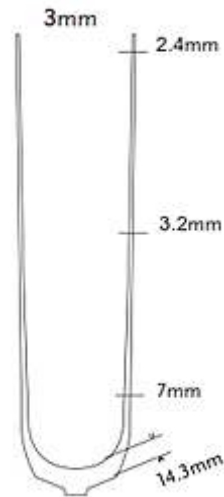

SIZING GUIDE

Measure the circumference of the stump at 4 cm above the distal end.

Measured sized (cm)	Locking 3mm Uniform
22-23.5 cm	FD1222
23.5-25 cm	FD1223.5
25-26.5 cm	FD1225
26.5-28 cm	FD1226
28-30 cm	FD1228
30-32 cm	FD1230
32-34 cm	FD1232
34-36 cm	FD1234
36-38 cm	FD1236
38-40 cm	FD1238

CristalSil Liner – Locking



The CristalSil liner is designed with silicone gel and without a fabric cover for interim prosthesis. The cristalSil liner is indicated for first Interim prosthesis after transtibial amputation.

Features and Benefits:

- Easy Glide surface coating for easy doning and doffing
- Medical grade silicone gel to protect the skin and offer the maximum comfort
- Easy to clean



Indications

Mobility level	Low to moderate activity
Amputation level	Transtibial


General information

Profile	3 mm
Suspension Method	Locking
Matrix	10 cm
Sizes	22,23.5,25,26.5,28,30,32,34,36,38,40,
Article Number	FD11xx

CristalSil Liner – Ordering Guidelines

PROFILE

3 mm uniform

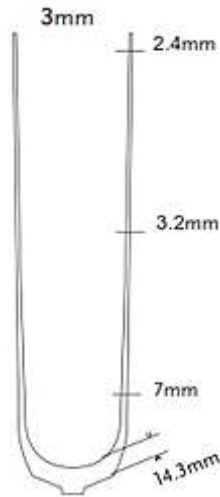

SIZING GUIDE

Measure the circumference of the stump at 4 cm above the distal end.

Measured sized (cm)	3mm Uniform
22-23.5 cm	FD11322
23.5-25 cm	FD11323.5
25-26.5 cm	FD11325
26.5-28 cm	FD11326
28-30 cm	FD11328
30-32 cm	FD11330
32-34 cm	FD11332
34-36 cm	FD11334
36-38 cm	FD11336
38-40 cm	FD11338

TruGel Liner – Locking



The TruGel liner is designed with very soft medical grade thermoplastic elastomer (TPE) and mineral oil. The improved gel formula provides to users durability and in the same time ensures maximum comfort. The one-way stretch fabric minimize pistoning and ensure flexibility and durability. The TruGel liner is indicated for users with Sensitive residual limb.

Features and Benefits:

- Soft gel to offer the best comfort
- Medical grade gel for care for the skin
- Durable fabric cover for easy doning and doffing
- Flexible distal cup
- Silky internal surface

Indications

Mobility level	Low to High activity
Amputation level	Transtibial

General information


Profile	6 mm – 3 mm
Suspension Method	Locking
Matrix	10 cm
Sizes	18,22,24,26,28,32,38,45
Article Number	• FDG146xx (6mm profile) • FDG143xx (3mm profile)




TruGel Liner – Ordering Guidelines

PROFILE

6 mm uniform

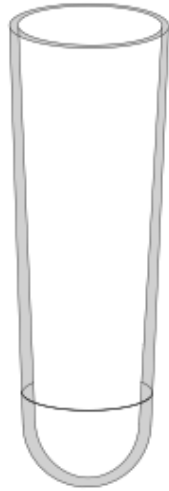

3 mm uniform

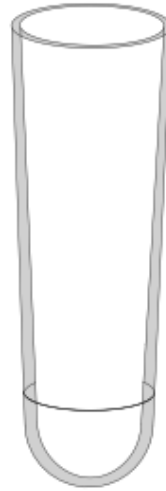

SIZING GUIDE

Measure the circumference of the stump at 4 cm above the distal end.

Measured sized (cm)	3mm Uniform	6mm Uniform
20-23 cm	FDG14322	FDG14620
24-25 cm	FDG14324	FDG14624
26-27 cm	FDG14326	FDG14626
28-31 cm	FDG14328	FDG14628
32-37 cm	FDG14332	FDG14632
38-43 cm	FDG14338	FDG14638
44-55 cm	FDG14344	FDG14644

TruGel Liner – Cushion



The TruGel liner is designed with very soft medical grade thermoplastic elastomer (TPE) and mineral oil. The improved gel formula provides to users durability and in the same time ensures maximum comfort. The one-way stretch fabric minimize pistoning and ensure flexibility and durability. The TruGel liner is indicated for users with Sensitive residual limb.

Features and Benefits:

- Soft gel to offer the best comfort
- Medical grade gel for care for the skin
- Durable fabric cover for easy doning and doffing
- Flexible distal cup
- Silky internal surface

Indications

Mobility level	Low to High activity
Amputation level	Transtibial

General information


Profile	6 mm
Suspension Method	Cushion
Matrix	without
Sizes	18,22,24,26,28,32,38,45
Article Number	• FDG046xx (6mm profile) • FDG043xx (3mm profile)




TruGel Liner – Ordering Guidelines

PROFILE

6 mm uniform

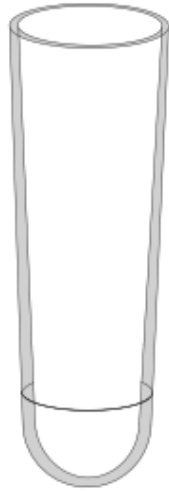

3 mm uniform

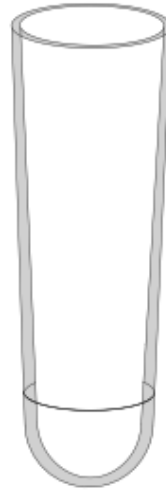

SIZING GUIDE

Measure the circumference of the stump at 4 cm above the distal end.

Measured sized (cm)	3mm Uniform	6mm Uniform
20-23 cm	FDG04322	FDG14620
24-25 cm	FDG04324	FDG14624
26-27 cm	FDG04326	FDG14626
28-31 cm	FDG04328	FDG14628
32-37 cm	FDG04332	FDG14632
38-43 cm	FDG04338	FDG14638
44-55 cm	FDG04344	FDG14644

TruGel Pro Liner – Locking



The TruGel Pro liner is designed with very soft medical grade thermoplastic elastomer (TPE) and mineral oil. The improved gel formula with skin care ingredients protect the skin and ensures maximum comfort. The one-way stretch fabric minimize pistoning and ensure flexibility and durability. The TruGel Pro liner is indicated for users with diabetic and allergic skin problem.

Features and Benefits:

- Soft gel to offer the best comfort
- skin care ingredients protect the skin
- Flexible distal cup
- Silky internal surface



TPE Gel Liners - Transtibial

Indications

Mobility level	Low to High activity
Amputation level	Transtibial

General information



Profile	6 mm – 3 mm
Suspension Method	Locking
Matrix	10 cm
Sizes	18,22,24,26,28,32,38,45
Article Number	• FDG136xx (6mm profile) • FDG133xx (3mm profile)

TruGel Pro Liner – Ordering Guidelines

PROFILE

6 mm uniform

3 mm uniform

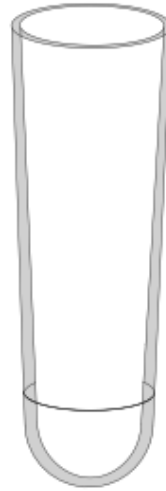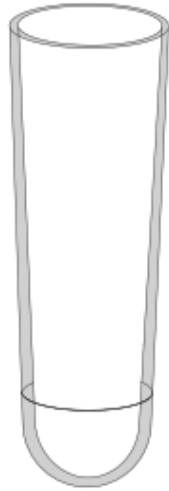

SIZING GUIDE

Measure the circumference of the stump at 4 cm above the distal end.

Measured sized (cm)	3mm Uniform	6mm Uniform
20-23 cm	FDG13322	FDG13620
24-25 cm	FDG13324	FDG13624
26-27 cm	FDG13326	FDG13626
28-31 cm	FDG13328	FDG13628
32-37 cm	FDG13332	FDG13632
38-43 cm	FDG13338	FDG13638
44-55 cm	FDG13344	FDG13644

TruGel pro Liner – Cushion



The TruGel Pro liner is designed with very soft medical grade thermoplastic elastomer (TPE) and mineral oil. The improved gel formula with skin care ingredients protect the skin and ensures maximum comfort. The one-way stretch fabric minimize pistoning and ensure flexibility and durability. The TruGel Pro liner is indicated for users with diabetic and allergic skin problem.

Features and Benefits:

- Soft gel to offer the best comfort
- skin care ingredients protect the skin
- Flexible distal cup
- Silky internal surface

Indications

Mobility level	Low to High activity
Amputation level	Transtibial

General information


Profile	6 mm – 3 mm
Suspension Method	Cushion
Matrix	without
Sizes	18,22,24,26,28,32,38,45
Article Number	• FDG036xx (6mm profile) • FDG033xx (3mm profile)




TruGel Pro Liner – Ordering Guidelines

PROFILE

6 mm uniform

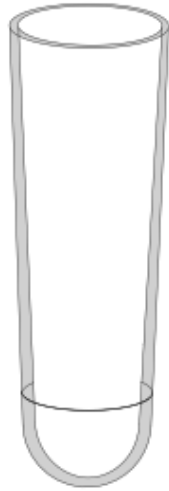

3 mm uniform

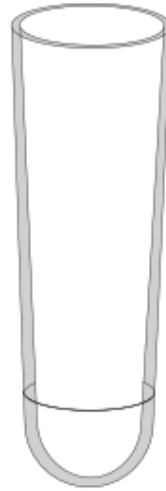

SIZING GUIDE

Measure the circumference of the stump at 4 cm above the distal end.

Measured sized (cm)	3mm Uniform	6mm Uniform
20-23 cm	FDG03322	FDG13620
24-25 cm	FDG03324	FDG13624
26-27 cm	FDG03326	FDG13626
28-31 cm	FDG03328	FDG13628
32-37 cm	FDG03332	FDG13632
38-43 cm	FDG03338	FDG13638
44-55 cm	FDG03344	FDG13644

GoSeal Liner – Cushion



The GoSeal liner is designed with very soft medical grade thermoplastic elastomer (TPE) and mineral oil. The GoSeal liner is indicated for standard use and offer Secur suspension by enhancing the sealing of the liner onto the sleeve.

Features and Benefits:

- Soft gel to offer the best comfort
- Medical grade gel for care for the skin
- Improve the sealing of the liner onto the sleeve
- Durable fabric cover for easy doning and doffing
- Silky internal surface



Indications

Mobility level	Low to Moderate activity
Amputation level	Transtibial

General information



Profile	6 mm - 3mm
Suspension Method	Cushion
Matrix	Without
Sizes	18,22,24,26,28,32,38,45
Article Number	• FDG006xx (6mm profile) • FDG003xx (3mm profile)

GoSeal Liner – Ordering Guidelines

PROFILE

6 mm uniform

3 mm uniform

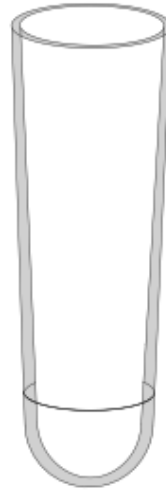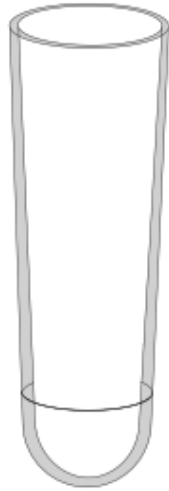

SIZING GUIDE

Measure the circumference of the stump at 4 cm above the distal end.

Measured sized (cm)	Locking 3mm Uniform	Locking 6mm Uniform
20-23 cm	FDG00322	FDG00620
24-25 cm	FDG00324	FDG00624
26-27 cm	FDG00326	FDG00626
28-31 cm	FDG00328	FDG00628
32-37 cm	FDG00332	FDG00632
38-43 cm	FDG00338	FDG00638
44-55 cm	FDG00344	FDG00644

EasySeal TF Silicone Liner



The EasySeal TF liner is designed for transfemoral sealing suction or elevated vacuum suspension. The EasySeal Liner is made with medical grade silicone offering comfort and durability to the user. The 4 suspension rings ensure secure suspension. The Easy Glide coating offers easy donning and easy clean.

Features and Benefits:

- Skin care formulation for care for the skin
- Fabric less exterior with special Easy Glide coating for easy donning and doffing
- Silky internal surface



Indications

Mobility level	Low to Active activity
Amputation level	Transfemoral


General information

Profile	3 mm
Suspension Method	Cushion
Matrix	Full matrix
Sizes	25,26.5,28,30, 32,34,36,38,40,45,50
Article Number	FDTF04XX

EasySeal TF Liner – Ordering Guidelines

PROFILE

3 mm uniform

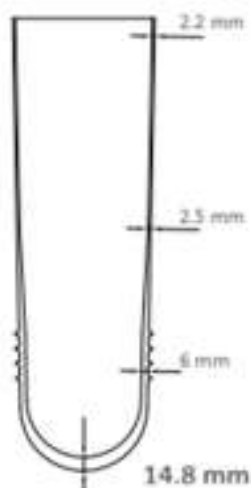

SIZING GUIDE

Measure the circumference of the stump at 4 cm above the distal end.

Measured sized (cm)	3mm Uniform
28-30 cm	FDTF0428
30-32 cm	FDTF0430
32-34 cm	FDTF0432
34-36 cm	FDTF0434
36-38 cm	FDTF0436
38-40 cm	FDTF0438
40-44 cm	FDTF0440
45-50 cm	FDTF0445

TF Silicone Liner – Locking



The TF liner is designed for transfemoral amputation. The TF liner offers to the user stability, comfort and security suspension. The TF liner includes many skin care ingredients to protect the skin and to ensure to the user the maximum comfort.

Features and Benefits:

- Skin care Silicone formulation for care for the skin
- Durable fabric cover for easy doning and doffing
- Flexible distal cup
- Silky internal surface

Indications

Mobility level	Low to moderate activity
Amputation level	Transfemoral

General information

Profile	3 mm
Suspension Method	Locking
Matrix	Full matrix
Sizes	28,30, 32,34,36,38,40,45,50
Article Number	FDTF1XX



TF Silicone Liner – Ordering Guidelines

PROFILE

3 mm uniform



SIZING GUIDE

Measure the circumference of the stump at 4 cm above the distal end.

Measured sized (cm)	3mm Uniform
28-30 cm	FDTF128
30-32 cm	FDTF130
32-34 cm	FDTF132
34-36 cm	FDTF134
36-38 cm	FDTF136
38-40 cm	FDTF138
40-44 cm	FDTF140
45-50 cm	FDTF145

TF Silicone Liner – Cushion



The TF liner is designed for transfemoral amputation. The TF liner offers to the user stability, comfort and security suspension. The TF liner includes many skin care ingredients to protect the skin and to ensure to the user the maximum comfort.

Features and Benefits:

- Skin care formulation for care for the skin
- Durable fabric cover for easy doning and doffing
- Flexible distal cup
- Silky internal surface

Indications

Mobility level	Low to moderate activity
Amputation level	Transfemoral

General information


Profile	3 mm
Suspension Method	Cusion
Matrix	Without
Sizes	28,30, 32,34,36,38,40,45,50
Article Number	FDTF0xx



TF Silicone Liner – Ordering Guidelines

PROFILE

3 mm uniform

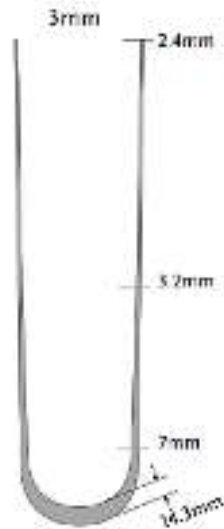

SIZING GUIDE

Measure the circumference of the stump at 4 cm above the distal end.

Measured sized (cm)	3mm Uniform
28-30 cm	FDTF028
30-32 cm	FDTF030
32-34 cm	FDTF032
34-36 cm	FDTF034
36-38 cm	FDTF036
38-40 cm	FDTF038
40-44 cm	FDTF040
45-50 cm	FDTF045

Basic Sleeve



The Basic sleeve is composed of TPE gel and a fabric cover. The Basic sleeve can be used in conjunction with any of our liners to ensure good security suspension.

Features and Benefits:

- Fabric cover and medical grade Gel
- Ensure perfect suspension
- Four Sizes



Indications

Mobility level	Low to Moderate activity
Amputation level	Transtibial and Transfemoral

General information

Profile	3 mm
Suspension Method	Cushion
Matrix	Without
Sizes	S, M, L, XL
Article Number	• FDS01xx

Basic Sleeve – Ordering Guidelines

PROFILE

3 mm uniform



SIZING GUIDE

Measure the circumference of the leg at 15 cm above the centre of the knee.

Proximal Circumference	Distal Circumference	3mm Uniform
28-36 cm	23-28 cm	FDS01S
36-46 cm	28-36 cm	FDS01M
43-56 cm	33-43 cm	FDS01L
51-69 cm	43-53 cm	FDS01XL

ProX Sleeve



The ProX sleeve is composed of TPE gel and a fabric cover. The ProFlex sleeve can be used in conjunction with any of our liners to ensure good security suspension. The ProX sleeve provides to the users comfort, stability and high performance for active users.

Features and Benefits:

- Fabric cover and medical grade Gel
- Ensure perfect suspension
- Seven Sizes



Knee Sleeve

Indications

Mobility level	Low to High activity
Amputation level	Transtibial and Transfemoral

General information


Profile	3 mm - 6 mm
Suspension Method	Cushion
Matrix	Without
Sizes	20,24,26,28,32,38,44
Article Number	• FDS023xx (3mm Profile) • FDS026xx (6mm Profile)

ProX Sleeve– Ordering Guidelines

PROFILE

3 mm uniform

6 mm uniform

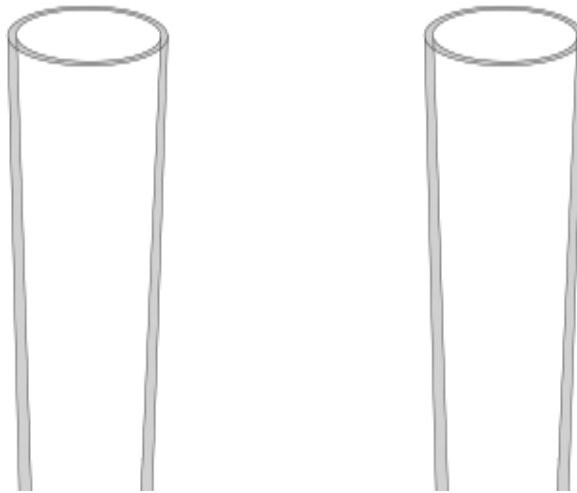

SIZING GUIDE

Measure the circumference of the leg at 15 cm above the centre of the knee.

Measured Size (cm)	3mm Uniform	6 mm Uniform
20-24 cm	FDS02320	FDS02620
25-30 cm	FDS02324	FDS02620
31-34 cm	FDS02326	FDS02620
35-37 cm	FDS02328	FDS02620
38-40 cm	FDS02332	FDS02620
41-45 cm	FDS02338	FDS02620
46-60 cm	FDS02326	FDS02620

Galaxy C4 LP Carbon foot



The Low profile Galaxy C4 LP is engineered to (K4) active users with limited clearance and provides durability and stability. The Galaxy C4 LP offers to the user both comfortable walking and dynamic rollover and offers walking at the highest level.

Features and Benefits:

- High energy return
- Progressive forefoot dynamics
- low profile design for limited clearance
- Lightweight and smooth gait
- Weatherfriendly
- Rated for 136 kg
- For low activity level



Galaxy C4 LP Carbon foot

Indications

Mobility level	Low to High activity
Amputation level	Transtibial and transfemoral
Max. body weight	166 kg

General information

Heel Height	10 mm
Connector	Male pyramid
Sizes	22 cm,23 cm,24 cm,25 cm,26 cm,27 cm,28 cm,29 cm,30 cm
Build Height (size 27)	68 mm (include pyramid and foot shell)
Foot Weight	560 g (include pyramid and foot shell)
Categories	Sizes 22-24; Cat. 1-6 Sizes 25-30; Cat. 1-9
Warranty	36 months for carbon foot 6 months for foot shell
Heel Wedge Hardness	Soft, Hard
Article Number	FD4CFLP



Galaxy C4 LP Carbon foot

Determining the Hardness Category

Weight kg	45-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147	148-166
Low Activity	1	1	2	3	4	5	6	7	8	9
Moderate Activity	1	2	3	4	5	6	7	8	9	N/A
High Activity	2	3	4	5	6	7	8	9	N/A	N/A

Galaxy C4 LP Carbon foot

Accessories

Foot Shell CS4

PU foot shell for Galaxy C4 LP foot with 6 months warranty with beige and brown color option



Foot Shell Informations

Sizes	22	23	24	25	26	27	28	29	30
Weight	150 g	155 g	195 g	210 g	226 g	235 g	265 g	275 g	310 g
Sides	Left (L), Right (R)								
Colors	Beige (C), Brown (B)								
Article N.	CS2								

Spectra Socks SS4

Protective spectra socks to protect the foot and to reduces the noise between foot and the foot shell



Spectra Socks informations

Article Number	SS4
Sizes	Small, Large
Color	Black

Heel Wedge Set

Article Number	Feet Size
SS40	22-27
SS41	28-30

Galaxy C4 LP Carbon foot

Ordering Guide

Foot Module Part Kit



Article N.	Size 22-30	Category 1-6	Side L-Left R-Right	Shell color C-Beige B-Brown

Foot Module Part Number



Article N.	Size 22-30	Category 1-6	Side L-Left R-Right

Foot Shell Part



Article N.	Size 22-30	Category 1-6	Side L-Left R-Right

Galaxy C3 Carbon foot



The Galaxy C3 is engineered to (K3) moderate users for many indoor and outdoor daily activities. The Galaxy C3 is designed to provide users with both comfortable walking and dynamic rollover. The light weight of the Galaxy C3 makes it the ideal foot for users seeking for stability and natural waking.

Features and Benefits:

- High energy return
- Progressive forefoot dynamics
- Lightweight and smooth gait
- Split toe for more stability and security
- Waterproof
- Weatherfriendly
- Rated for 136 kg
- For low to moderate activity level



Galaxy C3 Carbon foot

Indications

Mobility level	Moderate to High activity
Amputation level	Transtibial and transfemoral
Max. body weight	136 kg

General information

Heel Height	10 mm
Connector	Male pyramid
Sizes	22 cm,23 cm,24 cm,25 cm,26 cm,27 cm,28 cm,29 cm,30 cm
Build Height (size 27)	130mm (include pyramid and foot shell)
Foot Weight	520 g (include pyramid and foot shell)
Categories	Sizes 22-24; Cat. 1-4 Sizes 25-30; Cat. 1-6
Warranty	36 months for carbon foot 6 months for foot shell
Heel Wedge Hardness	Soft, Hard
Article Number	FD3CF



Galaxy C3 Carbon foot

Determining the Hardness Category

Weight kg	45-61	62-80	81-97	98-114	115-136
Low Activity	1	2	3	4	5
Moderate Activity	2	2	3	4	5

Galaxy C3 Carbon foot

Accessories

Foot Shell CS3

PU foot shell for Galaxy C3 foot with 6 months warranty with beige and brown color option



Foot Shell Informations

Sizes	22	23	24	25	26	27	28	29	30
Weight	150 g	155 g	195 g	210 g	226 g	235 g	265 g	275 g	310 g
Sides	Left (L), Right (R)								
Colors	Beige (C), Brown (B)								
Article N.	CS3								

Spectra Socks SS3

Protective spectra socks to protect the foot and to reduces the noise between foot and the foot shell



Spectra Socks informations

Article Number	SS3
Sizes	Small, Large
Color	Black

Heel Wedge Set

Article Number	Feet Size
SS40	22-27
SS41	28-30

Galaxy C3 Carbon foot

Ordering Guide

Carbon Feet

Foot Module Part Kit



Article
N.

Size
22-30

Category
1-6

Side
L-Left
R-Right

Shell color
C-Beige
B-Brown

Foot Module Part Number



Article
N.

Size
22-30

Category
1-6

Side
L-Left
R-Right

Foot Shell Part



Article
N.

Size
22-30

Category
1-6

Side
L-Left
R-Right

Galaxy C4 Symes Carbon foot



The Galaxy C3 Symes is engineered to (K4) high to moderate users with symes and pyrogoff amputations or transtibial users with long residual limbs for many indoor and outdoor daily activities. The Galaxy C3 Symes is designed to provide users with both comfortable walking and dynamic rollover. The light weight of the Galaxy C3 makes it the ideal foot for users seeking for stability and natural walking.

Features and Benefits:

- High energy return
- Progressive forefoot dynamics
- Lightweight and smooth gait
- Split toe for more stability and security
- Waterproof
- Weatherfriendly
- Rated for 166 kg
- For low to moderate activity level



Galaxy C4 Symes Carbon foot

Indications

Mobility level	High to Moderate to High activity
Amputation level	Transtibial and transfemoral
Max. body weight	166 kg

General information

Heel Height	10 mm
Connector	Lamination adapter
Sizes	22 cm,23 cm,24 cm,25 cm,26 cm,27 cm,28 cm,29 cm,30 cm
Build Height (size 27)	55 mm (includes foot shell)
Foot Weight	545 g (include pyramid and foot shell)
Categories	Sizes 22-24; Cat. 1-6 Sizes 25-30; Cat. 4-9
Warranty	36 months for carbon foot 6 months for foot shell
Heel Wedge Hardness	Soft, Hard
Article Number	FD3SCF



Galaxy C4 Symes Carbon foot

Determining the Hardness Category

Weight kg	45-52	53-59	60-68	69-77	78-88	89-100	101-116	117-130	131-147	148-166
Low Activity	1	1	2	3	4	5	6	7	8	9
Moderate Activity	1	2	3	4	5	6	7	8	9	N/A
High Activity	2	3	4	5	6	7	8	9	N/A	N/A

Galaxy C4 Symes Carbon foot

Accessories

Foot Shell CS4

PU foot shell for Galaxy C4 foot with 6 months warranty with beige and brown color option



Foot Shell Informations

Sizes	22	23	24	25	26	27	28	29	30
Weight	150 g	155 g	195 g	210 g	226 g	235 g	265 g	275 g	310 g
Sides	Left (L), Right (R)								
Colors	Beige (C), Brown (B)								
Article N.	CS3								

Spectra Socks SS4

Protective spectra socks to protect the foot and to reduces the noise between foot and the foot shell



Spectra Socks informations

Article Number	SS4
Sizes	Small, Large
Color	Black

Heel Wedge Set

Article Number	Feet Size
SS40	22-27
SS41	28-30

Galaxy C4 symes Carbon foot

Ordering Guide

Carbon Feet

Foot Module Part Kit



Article
N.

Size
22-30

Category
1-6

Side
L-Left
R-Right

Shell color
C-Beige
B-Brown

Foot Module Part Number



Article
N.

Size
22-30

Category
1-6

Side
L-Left
R-Right

Foot Shell Part



Article
N.

Size
22-30

Category
1-6

Side
L-Left
R-Right

Galaxy C2 Carbon foot



The Galaxy C2 is a lightweight foot engineered to (K2) low to moderate activity users for day to day activities. The Galaxy C2 is designed to provide users with stability and durability. The low profile design meets the needs for the users with limited clearance for indoor and outdoor environments

Features and Benefits:

- High energy return
- Progressive forefoot dynamics
- low profile design for limited clearance
- Lightweight and smooth gait
- Weatherfriendly
- Rated for 136 kg
- For low activity level



Galaxy C2 Carbon foot

Indications

Mobility level	K2 activity level
Amputation level	Transtibial and transfemoral
Max. body weight	136 kg

General information

Heel Height	10 mm
Connector	Male pyramid
Sizes	22 cm,23 cm,24 cm,25 cm,26 cm,27 cm,28 cm,29 cm,30 cm
Build Height (size 27)	68 mm (include pyramid and foot shell)
Foot Weight	490 g (include pyramid and foot shell)
Categories	1-5
Warranty	24 months for carbon foot 6 months for foot shell
Heel Wedge Hardness	Soft, Hard
Article Number	FD2CF



Galaxy C2 Carbon foot

Determining the Hardness Category

Weight kg	45-61	62-80	81-97	98-114	115-136
Low Activity	1	2	3	4	5
Moderate Activity	2	2	3	4	5

Galaxy C2 Carbon foot

Accessories

Foot Shell CS2

PU foot shell for Galaxy C2 foot with 6 months warranty with beige and brown color option



Foot Shell Informations

Sizes	22	23	24	25	26	27	28	29	30
Weight	150 g	155 g	195 g	210 g	226 g	235 g	265 g	275 g	310 g
Sides	Left (L), Right (R)								
Colors	Beige (C), Brown (B)								
Article N.	CS2								

Spectra Socks SS2

Protective spectra socks to protect the foot and to reduces the noise between foot and the foot shell



Spectra Socks informations

Article Number	SS2
Sizes	Small, Large
Color	Black

Heel Wedge Set

Article Number	Feet Size
SS40	22-27
SS41	28-30

Galaxy C2 Carbon foot

Ordering Guide

Foot Module Part Kit



Article
N.

Size
22-30

Category
1-6

Side
L-Left
R-Right

Shell color
C-Beige
B-Brown

Foot Module Part Number



Article
N.

Size
22-30

Category
1-6

Side
L-Left
R-Right

Foot Shell Part



Article
N.

Size
22-30

Category
1-6

Side
L-Left
R-Right

Galaxy C1 Carbon foot



The Galaxy C1 is a lightweight foot engineered to (K2) low to moderate activity users for daily activities. The foam heel provides soft and stable heel strike and smooth rollover.

Features and Benefits:

- High energy return
- Progressive forefoot dynamics
- Lightweight and smooth gait
- Weatherfriendly
- Rated for 136 kg
- For low activity level



Galaxy C1 Carbon foot

Indications

Mobility level	K2 activity level
Amputation level	Transtibial and transfemoral
Max. body weight	125 kg

General information

Heel Height	10 mm
Connector	Male pyramid
Sizes	22 cm,23 cm,24 cm,25 cm,26 cm,27 cm,28 cm,29 cm,30 cm
Build Height (size 27)	125mm (include pyramid and foot shell)
Foot Weight	570 g (include pyramid and foot shell)
Categories	1-5
Warranty	24 months for carbon foot 6 months for foot shell
Heel Wedge Hardness	Soft, Hard
Article Number	FD1CF



Carbon Foot

Galaxy C1 Carbon foot

Determining the Hardness Category

Weight kg	45-61	62-80	81-97	98-114	115-136
Low Activity	1	2	3	4	5
Moderate Activity	2	2	3	4	5

SACH Foot Geriatric

Indications

Mobility level	K2 activity level
Amputation level	Transtibial and transfemoral
Max. body weight	100 kg

General information

Heel Height	10 mm
Sizes	22 cm,23 cm,24 cm,25 cm,26 cm,27 cm,28 cm,29 cm,30 cm
Build Height (size 27)	80 mm (include pyramid and foot shell)
Foot Weight	330 g (Size 26)
Warranty	12 months
Article Number	FD0SFG



Carbon Feet

Single Axis FootFoot

Indications

Mobility level	K2 activity level
Amputation level	Transtibial and transfemoral
Max. body weight	100 kg

General information

Heel Height	10 mm
Sizes	22 cm,23 cm,24 cm,25 cm,26 cm,27 cm,28 cm,29 cm,30 cm
Build Height (size 27)	50 mm (include pyramid and foot shell)
Foot Weight	350 g (Size 26)
Warranty	12 months
Article Number	FD1SF



Shuttle Lock

General information

System weight	Aprox. 280 g
Max. patient weight	125 kg
Activity level	1 to 4
Amputation Level	Transtibial
Article Number	2F01=10



Locking System

Clamp Adapter Ø 30 mm

General information

System weight	Aprox. 80 g
Max. patient weight	100 kg
Tube Ø	30 mm
Material	Stainless Steel
Article Number	2F12S



General information

System weight	Aprox. 70 g
Max. patient weight	100 kg
Tube Ø	30 mm
Material	Aluminum
Article Number	2F12A



Tube Clamp Adapter Ø 30 mm

General information

Article Number	2F13S=450	2F13S=220
System weight	315 g	198 g
Max. patient weight	100 kg	
Tube Ø	30 mm	
Material	Stainless Steel	



General information

Article Number	2F13A=450	2F13A=220
System weight	250 g	160 g
Max. patient weight	100 kg	
Tube Ø	30 mm	
Material	Aluminum	



4 Prongs Socket Adapter

General information

System weight	Aprox. 90 g
Max. patient weight	100 kg
Material	Stainless Steel
Article Number	2F13S



General information

System weight	Aprox. 70 g
Max. patient weight	100 kg
Material	Aluminum
Article Number	2F13A



Lamination Anchor with Pyramid

General information

System weight	Aprox. 130 g
Max. patient weight	150 kg
Material	Stainless Steel
Article Number	2F19



General information

System weight	Aprox. 92 g
Max. patient weight	125kg
Material	Stainless steel
Article Number	2F20



Lamination Anchor 3 Prongs rotatable

General information

System weight	Aprox. 170 g
Max. patient weight	125 kg
Material	Stainless Steel
Article Number	2F22



General information

System weight	Aprox. 180 g
Max. patient weight	125kg
Material	Stainless steel
Article Number	2F21





CREATIVE ORTHOPEADICS SRL

Via F. Ugoni 36, 25126 (BS) Italy

TEL. +393500577531

EMAIL. Info@creativeorthopedics.it

www.creativeorthopedics.it